


Dental trolleys

Prophy 700



Jelly

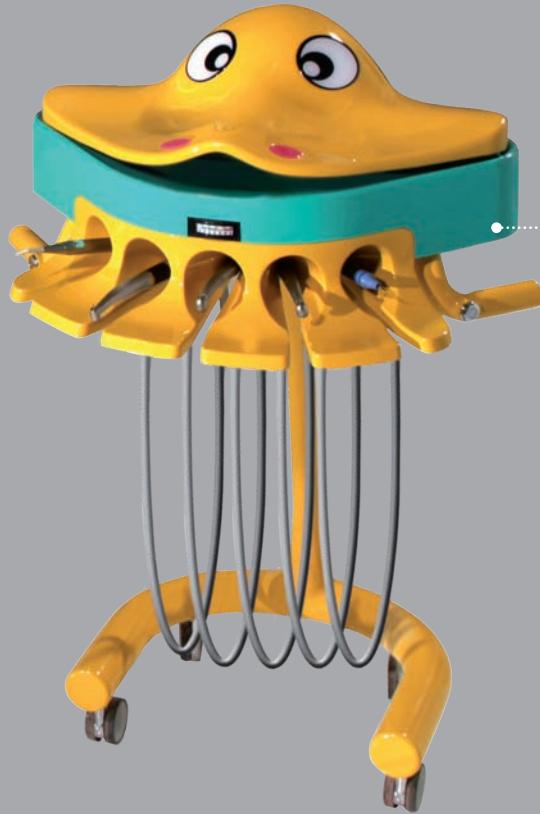

The Prophy 700 dental trolley is elegantly and ergonomically designed and permits total freedom of movement. The hanging cord instrument table can be customised to hold up to five electric or air driven instruments of your choice. Extremely lightweight and easy to move, the Prophy 700 promotes better coordinated interactions between dentist and assistant. The Jelly is a unique dental trolley dedicated to the wellbeing of your younger patients. Psychological research has shown that a young child's initial perceptions play a key role in determining his subsequent behaviour. Young patients are likely to perceive this trolley, with its funny jellyfish shape, as a toy rather than a threat. They will therefore see their visit to the dentist as a game, will not be afraid, and will see you as a friend to be trusted! Because they are fully independent operating units, these trolleys also allow you to carry on working if your main dental unit should suffer a breakdown.



All the main functions are within easy reach of your foot thanks to the semi-progressive foot control provided. The foot unit incorporates controls for instrument air drive and spray, micromotor reversal and endodontic scaler.



As an optional, the semi-progressive foot control that comes as standard equipment can be replaced by a Feed Fec foot control. The instrument foot control offers progressive action. It also lets you reverse the direction of micromotors and control scaler endodontic functions. Instrument spray can be switched on permanently or controlled by means of the lever at the bottom. Product code DE3.12FM



A satin finish polycarbonate base can be fitted as an option to both models of trolley. The base provides a convenient shelf for the foot control when the trolley has to be moved. Product code 0550325



The head of the jellyfish houses a built-in negatoscope as well as a generous 480x570mm tray.



An ABS instrument tray is also available for mounting over the equipment table of the Prophy 700. The tray is 375x275mm in size. Product code DE2.00M



All that is needed to make these trolleys fully functional is a compressed air supply and a low voltage socket – two items found on all dental units.

Data sheet and equipment list

	Jelly Product code DE800KP	Prophy 700 Product code DE700E	Product code of optional
Air supply	6/8 bar	6/8 bar	
Power supply	24/220 V - 50 Hz	24/220 V - 50 Hz	
Dimensions	60 x 100 x 55 cm	50 x 50 x 90 cm	
Weight	22.5 kg	15 kg	
Luzzani® Minimate steel syringe	●	●	
High speed handpiece cord	●3	●3	
Electric micromotor kit with optic fibre	Optional	Optional	Product code PM504FM
Ems® compatible piezoelectric scaler, 5 th instrument	Optional	Optional	Product code AG046BM
Satelec® compatible piezoelectric scaler, 5 th instrument	Optional	Optional	Product code AG055SM
Piezo One LED piezoelectric scaler, 5 th instrument		Optional	Product code AG047FM
LED curing light	Optional	Optional	Product code DEL1.01WM
Semi-progressive foot control	●	●	
Feed Fec progressive foot control	Optional	Optional	Product code DE3.12FM
100 cl spray water bottle	●	●	
Mains/bottle water switching kit	Optional	Optional	Product code PCI866AM
Satin finish polycarbonate base	Optional	Optional	Product code 0550325
Magnetic tray table		Optional	Product code DE2.00M
Negatoscope	●		
Warranty	24 months	24 months	